



Height adjustable over-table x-ray shield


model: 326/05



A This mobile shield combines all the attributes of overhead-suspended shields & table mounted shields in a single unit, without the need for installation. The front protective curtain can slide and is shown fully extended to the tableside.

B The counterbalanced window is shown at its maximum overall height of 177cm. When lowered it is 145cm high.

C The window corner cut-out can be positioned across the table and over the patient. The shield is easily adjusted for height & width, and provides full body protection without impeding procedures.

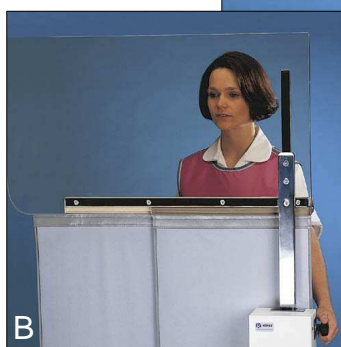

A dramatic reduction in radiation exposure to the whole body can be achieved when this mobile x-ray shield is properly utilised. It creates an extensive safety zone within which the practitioner can work. The wide shadow cast by the shield also provides protection to assisting personnel.

- The counterbalanced lead acrylic window may be positioned over the patient and height adjusted to suit the table.
- The 50x116cm window is supported at one end by a stand, which is shaped to allow for below table movements of C-arm mounted x-ray tubes.
- The 80cm high flexible curtains, which hang from the window, are adjustable for width from 45cm to 80cm. When extended to the tableside they provide maximum protection from scatter radiation.
- Lead equivalent: 0.5mm throughout.



Weights, dimensions and ordering information

Ordering information

Description: Height adjustable over-table x-ray shield.
Model: 326/05

Net weight: 51 kg

Shipping details

Quantity: 1 x heavy duty cardboard carton.
Gross weight: 60 kg
Packed size: 128 x 58 x 32cm

Weights & dimensions are approximate.

